




45 Gallon, 2 Shelves, 2 Door, Manual Close, Flammable Waste Cabinet, Sure-Grip® EX, White - 8945053

Product Images



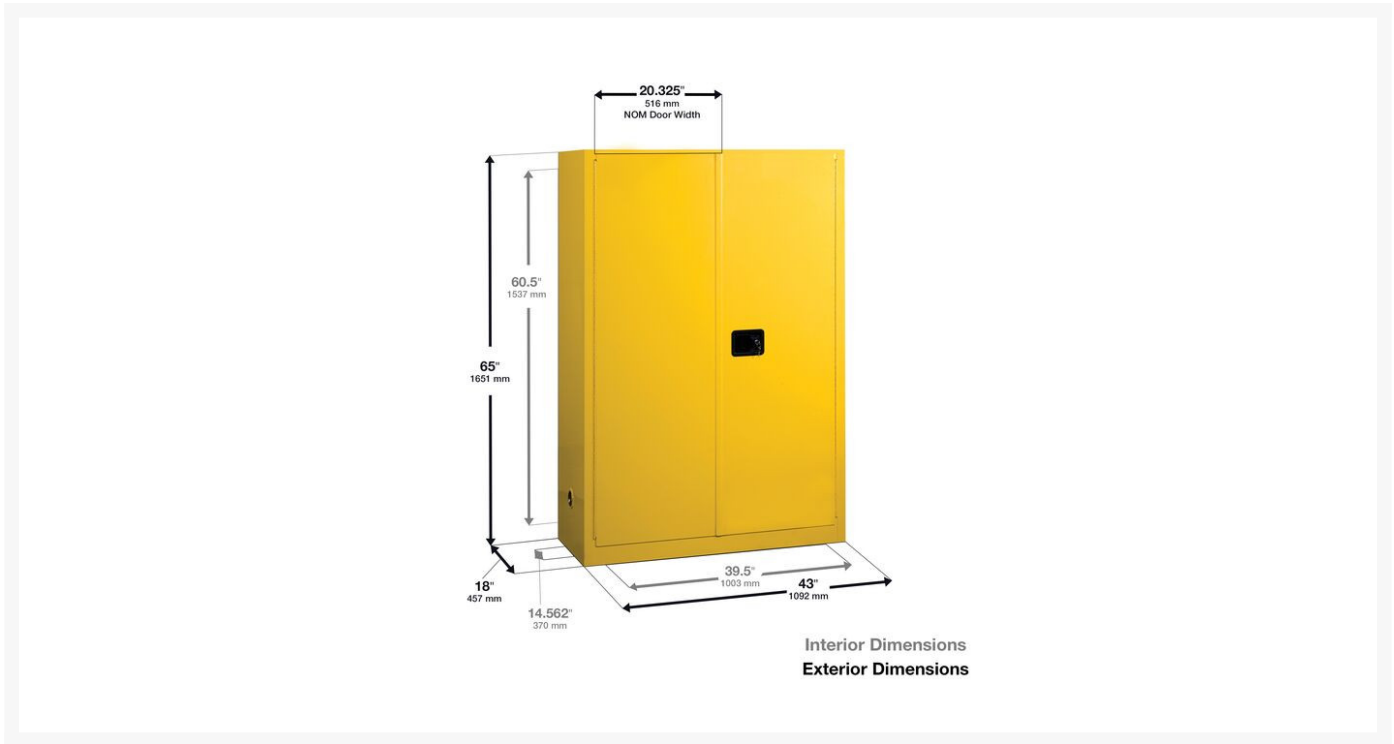

Short Description

- Ideal solution for temporary collection and segregation of waste liquids
- Large warning labels clearly identify contents
- FM-approved
- Ten-year warranty

Description

White Safety Cabinets for Flammable Waste are the solution for temporary collection and segregation of waste liquids. Large warning label "Flammable Waste" clearly identifies contents. Made to the same exacting standards as flammable safety cabinets, including double-wall construction, dual vents, grounding wire connections, adjustable shelves, leakproof sills, three-point self-latching doors and leveling feet.

Manual-close doors open to a full 180 degrees and self-latch when pushed closed. Cabinets are FM approved and come with a ten-year warranty.

Specifications

Model No	8945053
UPC	697841113794
International Model No.	89450531
Color	White
Pallet Quantity	2
Material Specifications	18-Gauge Steel
FM Approval	Yes
OSHA Compliance	Yes
Complies with NFPA Code 30	Yes
Capacity	45 gallon
Liter Capacity	170
Number of Doors	2
Door Type	Manual
Number of Shelves	2
Adjustable Shelves	Yes
Dimensions, Exterior	65" H x 43" W x 18" D
Dimensions, Exterior	1651mm H x 1092mm W x 457mm D
Dimensions, Interior	60.5" H x 39.5" W x 14.6" D
Dimensions, Interior	1537mm H x 1003mm W x 370mm D
Height (Inches)	65
Width (Inches)	43
Depth	18
Net Weight, lbs	270.0
Net Weight, kg	122.5