




# 45 Gallon, 2 Shelf, 2 Door, Self-Close, Hazardous Material Cabinet, ChemCor®, Royal Blue - 8645282

## Product Images



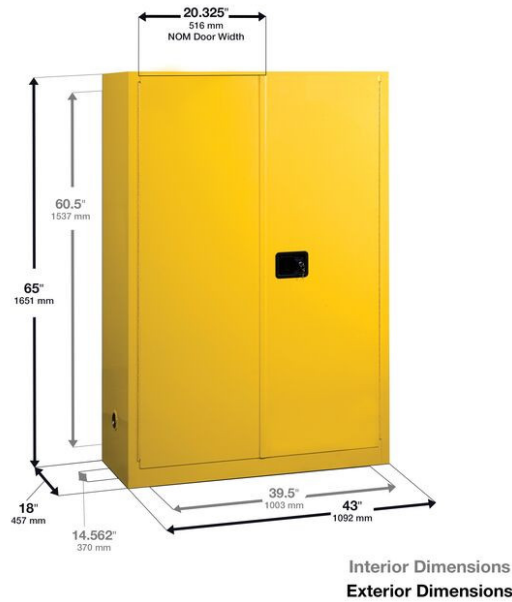

## Short Description

- Provide complete protection from harsh acids and corrosive chemicals
- Powder-coat finish provides increased chemical resistance
- Self-closing doors
- Seam-free surface eliminates the need for rivets or fasteners
- Complies with NFPA 1, NFPA 400 and the International Fire Code
- FM-approved
- Custom safety cabinets available, [get a quote here](#)

## Description

ChemCor® thermoplastic coating on all interior surfaces provides superior protection against corrosion from drips and off-gassing from harsh chemicals. Seam-free surface eliminates the need for rivets or fasteners which risk corrosion from seep-through chemicals. Non-porous surface resists staining and cleans up easily.

Looking for a custom safety cabinet? [Request a quote here](#) for your storage requirements.

ChemCor® Lined Hazardous Material Safety Cabinets have the same quality features found in Justrite's flammable cabinets — including double-wall, 18-gauge (1mm) welded steel construction, 2" (51mm) liquid-tight containment sump, dual vents, adjustable shelves, three-point stainless steel bullet self-latching doors with U•Loc™ padlockable handle and Haz-Alert™ reflective warning labels. An all epoxy baked-on powder-coat finish, inside and out, provides increased chemical resistance. Interior ChemCor® liner ensures complete corrosion resistance and protection. Self-close door(s) shuts and latches automatically when a fusible link melts at 165°F (74°C) under fire conditions. Unique, concealed self-closing mechanism offers obstruction-free

access to contents.

Exclusive to Justrite, each cabinet comes with a "Hazardous - Keep Fire Away" label, as well as a pack of ten application-specific labels to choose from and apply to the cabinet in order to indicate its contents.


Cabinets comply with NFPA 1, NFPA 400 Hazardous Material Code, and the International Fire Code. Cabinets are FM approved and come with a ten-year warranty.

# Specifications

Model No	8645282
UPC	697841154704
International Model No.	86452821
Color	Royal Blue
Pallet Quantity	2
Material Specifications	18-Gauge Steel
FM Approval	Yes
Complies with International Fire Code	Yes
Complies with NFPA Code 30	Yes
Complies with NFPA 1 Fire-Code	Yes
Complies with NFPA 400	Yes
Capacity	45 gallon
Liter Capacity	170
Number of Doors	2
Door Type	Self-close
Number of Shelves	2
Adjustable Shelves	Yes
Dimensions, Exterior	65" H x 43" W x 18" D
	1651mm H x 1092mm W x 457mm D
Dimensions, Interior	60.5" H x 39.5" W x 14.6" D
	1537mm H x 1003mm W x 370mm D
Height (Inches)	65
Width (Inches)	43
Depth	18
Net Weight, lbs	274.0
Net Weight, kg	124.3

## Prop 65

---

 **WARNING:** This product can expose you to chemicals including lead, which is known to the State of California to cause cancer and birth defects or other reproductive harm. For more information go to [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov).