




# 30 Gallon, 1 Shelf, 2 Doors, Manual Close, Corrosives/Acid Steel Safety Cabinet, Sure-Grip® EX, Blue - 893002

## Product Images



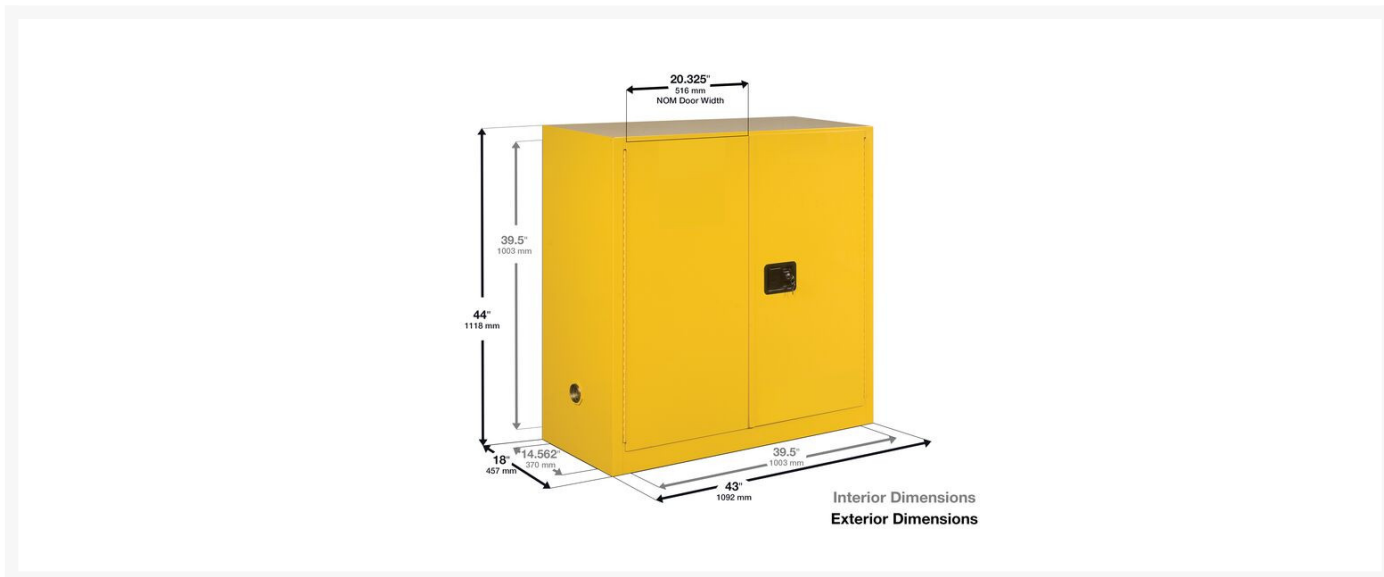

## Short Description

---

- Exclusively formulated paint and poly trays resist acids, bases and solvents
- Offer same quality safety features as cabinets for flammables
- FM-approved
- Ten-year warranty

## Description

---

Research environments using many different types of chemicals have special storage requirements. Whether using acids, bases or solvents, it is critical that incompatible liquids are segregated to avoid disastrous explosions or fire. The SDS (Safety Data Sheet) should be reviewed to determine chemical characteristics and recommended storage practices.

Justrite blue steel safety cabinets for corrosives have the same quality safety features as the cabinets for flammables. All models come with double-wall construction, dual vents, grounding wire connections, adjustable shelves, leakproof sills, three-point self-latching doors and leveling feet. Manual-close door(s) opens to a full 180 degrees and self-latches when pushed closed.

Additionally, they include polyethylene trays attached to galvanized steel shelves and a separate polyethylene liner for the bottom sump resist aggressive chemicals. Easily remove the liner for quick cleaning of drips and leaks. This popular 30-gallon (114-L) cabinet also includes an extra polyethylene work tray that can be secured to cabinet top for a handy work surface. An all epoxy baked-on powder-coat finish, inside and out, provides increased chemical resistance.

Manual close models meet OSHA and NFPA 30. Cabinets are FM approved and come with a ten-year warranty.

## Specifications

---

Model No	893002
UPC	697841113299
International Model No.	8930021
Color	Blue
Pallet Quantity	2
Material Specifications	18-Gauge Steel
FM Approval	Yes
OSHA Compliance	Yes
Complies with NFPA Code 30	Yes
Capacity	30 gallon
Liter Capacity	114
Number of Doors	2
Door Type	Manual
Number of Shelves	1
Adjustable Shelves	Yes
Dimensions, Exterior	44" H x 43" W x 18" D
	1118mm H x 1092mm W x 457mm D
Dimensions, Interior	39.6" H x 39.5" W x 14.6" D
	1006mm H x 1003mm W x 370mm D
Height (Inches)	44
Width (Inches)	43
Depth	18
Net Weight, lbs	196.5
Net Weight, kg	89.1