




90 Gallon, 2 Shelves, 2 Doors, Self Close, Corrosives/Acid Steel Safety Cabinet, Sure-Grip® EX, Blue - 899022

Product Images

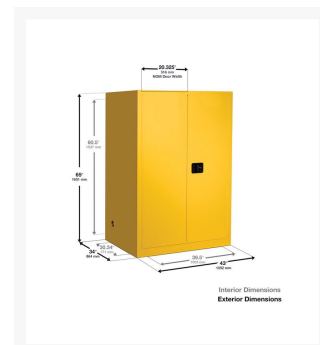

Short Description

- Exclusively formulated paint and poly trays resist acids, bases and solvents
- Offer same quality safety features as cabinets for flammables
- FM-approved
- Ten-year warranty

Description

Research environments using many different types of chemicals have special storage requirements. Whether using acids, bases or solvents, it is critical that incompatible liquids are segregated to avoid disastrous explosions or fire. The SDS (Safety Data Sheet) should be reviewed to determine chemical characteristics and recommended storage practices.

Justrite blue steel safety cabinets for corrosives have the same quality safety features as the cabinets for flammables. All models come with double-wall construction, dual vents, grounding wire connections, adjustable shelves, leakproof sills, three-point self-latching doors and leveling feet. Self-close door(s) shuts and latches automatically when a fusible link melts at 165°F (74°C) under fire conditions. Unique, concealed self-closing mechanism offers obstruction-free access to contents.

Additionally, they include polyethylene trays attached to galvanized steel shelves and a separate polyethylene liner for the bottom sump resist aggressive chemicals. Easily remove the liner for quick cleaning of drips and leaks. An all epoxy baked-on powder-coat finish, inside and out, provides increased chemical resistance.

Self Close models meet OSHA, NFPA 30, NFPA 1 and the International Fire Code. Cabinets are FM approved and come with a ten-year warranty.

Specifications

Model No	899022
UPC	697841114876
International Model No.	8990221
Color	Blue
Pallet Quantity	1
Material Specifications	18-Gauge Steel
FM Approval	Yes
OSHA Compliance	Yes
Complies with International Fire Code	Yes
Complies with NFPA Code 30	Yes
Complies with NFPA 1 Fire-Code	Yes
Capacity	90 gallon
Liter Capacity	341
Number of Doors	2
Door Type	Self-close
Number of Shelves	2
Adjustable Shelves	Yes
Dimensions, Exterior	65" H x 43" W x 34" D
	1651mm H x 1092mm W x 864mm D
Dimensions, Interior	60.5" H x 39.5" W x 30.3" D
	1537mm H x 1003mm W x 771mm D
Height (Inches)	65
Width (Inches)	43
Depth	34
Net Weight, lbs	387.0
Net Weight, kg	175.5

Prop 65

⚠️WARNING: This product can expose you to chemicals including lead, which is known to the State of California to cause cancer and birth defects or other reproductive harm. For more information go to www.P65Warnings.ca.gov.