




# 45 Gallon, 2 Shelves, 2 Doors, Manual Close, Corrosives/Acids Safety Cabinet, ChemCor®, Blue - 8945022

## Product Images



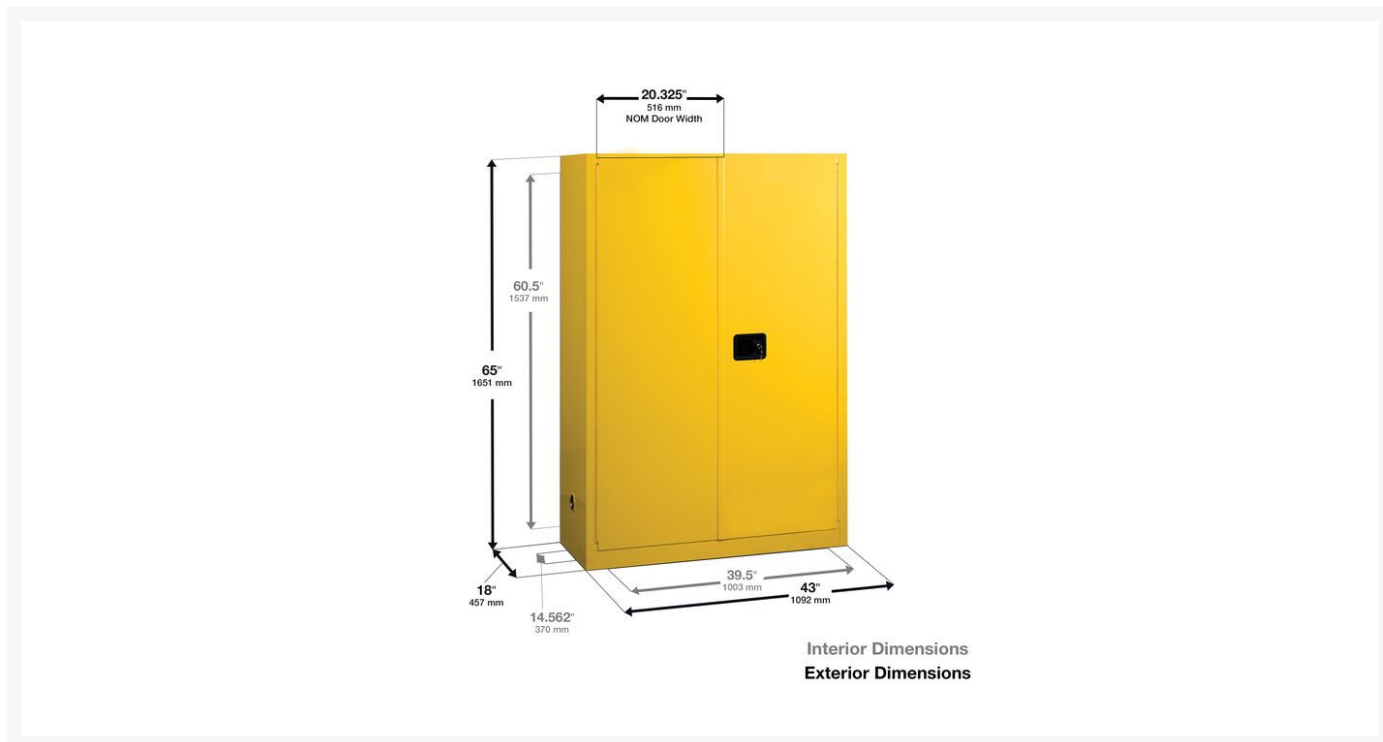

## Short Description

---

- Provide complete protection from harsh acids and corrosive chemicals
- Fire-protection and code compliance
- Seam-free surface eliminates the need for rivets or fasteners
- FM-approved

## Description

---

Research environments using many different types of chemicals have special storage requirements. Whether using acids, bases or solvents, it is critical that incompatible liquids are segregated to avoid disastrous explosions or fire. The SDS (Safety Data Sheet) should be reviewed to determine chemical characteristics and recommended storage practices.

The ChemCor® thermoplastic coating on all interior surfaces provides superior protection against corrosion from drips and off-gassing from harsh chemicals. Non-porous surface resists staining and cleans up easily.

ChemCor® Lined Acid Safety Cabinets have the same quality features found in Justrite's flammable cabinets — including double-wall, 18-gauge (1mm) welded steel construction, 2" (51mm) liquid-tight containment sump, dual vents, adjustable shelves, three-point stainless steel bullet self-latching doors with U•Loc™ padlockable handle and Haz-Alert™ reflective warning labels. An all epoxy baked-on powder-coat finish, inside and out, provides increased chemical resistance. Interior ChemCor® liner ensures complete corrosion resistance and protection. Manual-close door(s) opens to a full 180 degrees and self-latches when pushed closed.

Manual close models meet OSHA and NFPA 30. Cabinets are FM approved and come with a ten-year

warranty.

## Specifications

---

Model No	8945022
UPC	697841113749
International Model No.	89450221
Color	Blue
Pallet Quantity	2
Material Specifications	18-Gauge Steel
FM Approval	Yes
OSHA Compliance	Yes
Complies with NFPA Code 30	Yes
Capacity	45 gallon
Liter Capacity	170
Number of Doors	2
Door Type	Manual
Number of Shelves	2
Adjustable Shelves	Yes
Dimensions, Exterior	65" H x 43" W x 18" D
	1651mm H x 1092mm W x 457mm D
Dimensions, Interior	60.5" H x 39.5" W x 14.6" D
	1537mm H x 1003mm W x 370mm D
Height (Inches)	65
Width (Inches)	43
Depth	18
Net Weight, lbs	193.0
Net Weight, kg	87.5