


15 Gallon, 1 Shelf, 1 Door, Manual Close, Right Hinge, Corrosives/Acids Safety Cabinet, ChemCor® Under Fume Hood, Blue - 8825022

Product Images



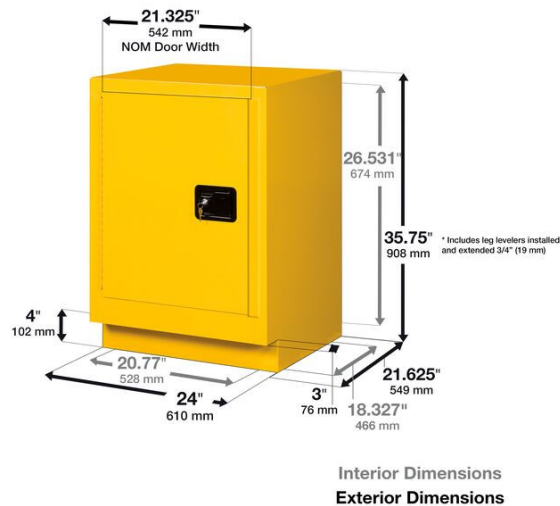

Short Description

- Provide complete protection from harsh acids and corrosive chemicals
- Fire-protection and code compliance
- Offer one adjustable shelf and kickplate for undercounter or fume hood base use
- Meet OSHA and NFPA regulations
- FM-approved
- Multiple color options

Description

Research environments using many different types of chemicals have special storage requirements. Whether using acids, bases or solvents, it is critical that incompatible liquids are segregated to avoid disastrous explosions or fire. The SDS (Safety Data Sheet) should be reviewed to determine chemical characteristics and recommended storage practices.

The ChemCor® thermoplastic coating on all interior surfaces provides superior protection against corrosion from drips and off-gassing from harsh chemicals. Seam-free surface eliminates the need for rivets or fasteners which risk corrosion from seep-through chemicals. Non-porous surface resists staining and cleans up easily.

ChemCor® Lined Under Fume Hood Cabinets has one adjustable shelf and a 4"h x 3"d (102mm x 77mm) toe kick for undercounter or fume hood base use. Single-door cabinets are offered in either left or right-hinge door closures for easy operation where space is at a premium. Back wall vents with built-in flame arresters are positioned at lower left and top right for jurisdictions that mandate venting of vapors. Specially-formulated, baked-on powder coat paint on cabinet exterior provides for long life and extra resistance to harsh chemicals. In addition, they have all the great features and benefits of blue steel cabinets, including grounding wire connections, adjustable shelves, leakproof sills, three-point self-latching doors and leveling

feet. Manual-close door opens to a full 180 degrees and self-latches when pushed closed.

Quality-backed by Justrite's Ten-Year Warranty. Cabinets meet OSHA and NFPA regulations, and are FM approved.

Specifications

Model No	8825022
UPC	697841120570
International Model No.	8825022
Color	Blue
Material Specifications	18-Gauge Steel
FM Approval	Yes
OSHA Compliance	Yes
Complies with NFPA Code 30	Yes
Capacity	15 gallon
Liter Capacity	57
Number of Doors	1
Door Type	Manual
Number of Shelves	1
Adjustable Shelves	Yes
Dimensions, Exterior	35.75" H x 24" W x 21.625" D
	908mm H x 610mm W x 549mm D
Dimensions, Interior	26.5" H x 20.8" W x 18.3" D
	674mm H x 528mm W x 466mm D
Height (Inches)	35.75
Width (Inches)	24
Depth	21.625
Net Weight, lbs	115.0
Net Weight, kg	52.2