


2.5 Gallon Steel Safety Can for Laboratories, Type I, Rigid Bottom Brass Faucet, Red - 7225150

Product Images



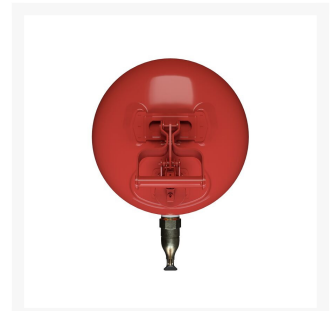

Short Description

- Ideal for laboratory settings
- Self-close, pressure relief caps
- Vacuum-relief valve within the lid
- Faucet mounted 2 inches above the bottom
- FM Approved

Description

Justrite's Steel Laboratory Safety Cans have all the benefits of the Type I cans. Added features make these the choice of lab professionals for safely dispensing hazardous liquids into test tubes, small flasks and beakers. Fill spouts have self-close, pressure-relief caps and stainless steel flame arresters. These cans also feature a vacuum-relief valve within the lid for smooth liquid flow, and self-closing faucets for trouble-free dispensing. Cans with the smaller No. 08540 faucet offer easier dispensing in lab environments, while the No. 08902 faucet is better suited for heavy, industrial use.

Shelf Cans have a faucet mounted 2 inches above the bottom to eliminate lifting the cans during dispensing.

FM Approved

Model No	7225150
UPC	697841140561
Brand	Justrite
International Model No.	7225150Z
FM Approval	Yes
OSHA Compliance	Yes
Complies with NFPA Code 30	Yes
Meets California Requirements	Required in California, CCR Title 8, § 1930(a)(1) CAL/OSHA; Exempt under CCR Title 13 § 2467.3(c)
Color	Red
Material Specifications	Steel/Brass Faucet
Dimensions, Exterior	10.6" H x 11.4" Outer Dia.
Dimensions, Exterior	270mm H x 290mm Outer Dia.
Capacity	2.5 gallon
Liter Capacity	9
Fill/Pour Hole inch	1.68
Fill/Pour Hole mm	43
Net Weight, lbs	5.9
Net Weight, kg	2.7

Specifications