
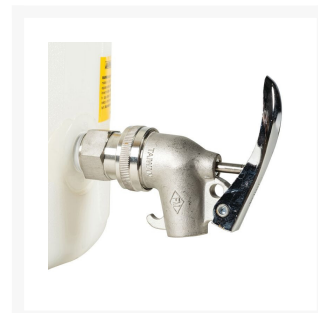


# 5 Gallon, Polypropylene Quick-Disconnect Disposal Safety Can with Faucet, Stainless Steel Fittings for 1/4" Tubing, White - 12770

## Product Images

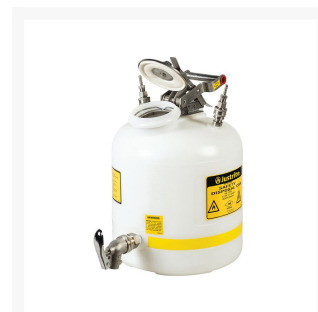

## Short Description

- Ideal for HPLC waste applications
- Provides hands-free solvent collection
- Durable and safe
- Features quick-disconnect system

## Description

---

Ideal for HPLC waste applications:

- Vapor and fire protection
- Safeguard against accidental spills
- Maximize safety and convenience

Hands-free solvent collection: Solvent flows freely from your process into the container through tubing making it suitable for either pump-fed or gravity-flow HPLC waste methods. Two disconnect fittings, one intake and one vent, eliminate back pressure in liquid flow disposal line and provide controlled handling of vapor.

Durability and safety: Translucent polyethylene container allows quick visual inspection of liquid level. Pouring spout is equipped with a stainless steel flame arrester, which reliably dissipates heat to prevent flashback ignition. Spring-loaded, sealed lid automatically vents to guard against explosion. Control odors with a coalescing/carbon filter.

### *Quick-Disconnect system*

- Easy-on, easy-off mechanism offers detachment of the can from the process without leaks, drips or vapor release. Once disconnected, the container is automatically sealed to prevent spills.
- Valved quick-disconnect fittings automatically close when decoupled to prevent vapor or liquid release.
- Stainless-steel fittings offer PTFE O-rings on standard containers. Also available with EPDM O-ring seals on customized cans.
- Coalescing carbon filterOptional filters are for use with only round style cans.
- Easy-carrying handle with fixed position latch keeps spring loaded cover open for convenient pouring during disposal.
- Economically priced polypropylene barbed fittings offer EPDM O-rings.
- Translucent color lets you check liquid level at a glance.

## Specifications

---

Model No	12770
UPC	697841102316
International Model No.	12770
FM Approval	Yes
OSHA Compliance	Yes
Complies with NFPA Code 30	Yes
Color	White
Material Specifications	Polyethylene
Dimensions, Exterior	20" H x 12" Outer Dia. 508mm H x 305mm Outer Dia.
Capacity	5 gallon
Liter Capacity	19
Fill/Pour Hole inch	3.65
Fill/Pour Hole mm	93
Net Weight, lbs	11.0
Net Weight, kg	5.0