


Wall Bracket Kit With Stainless Connector to Relocate The HPLC Coalescing/Carbon Filter - 28151

Product Images

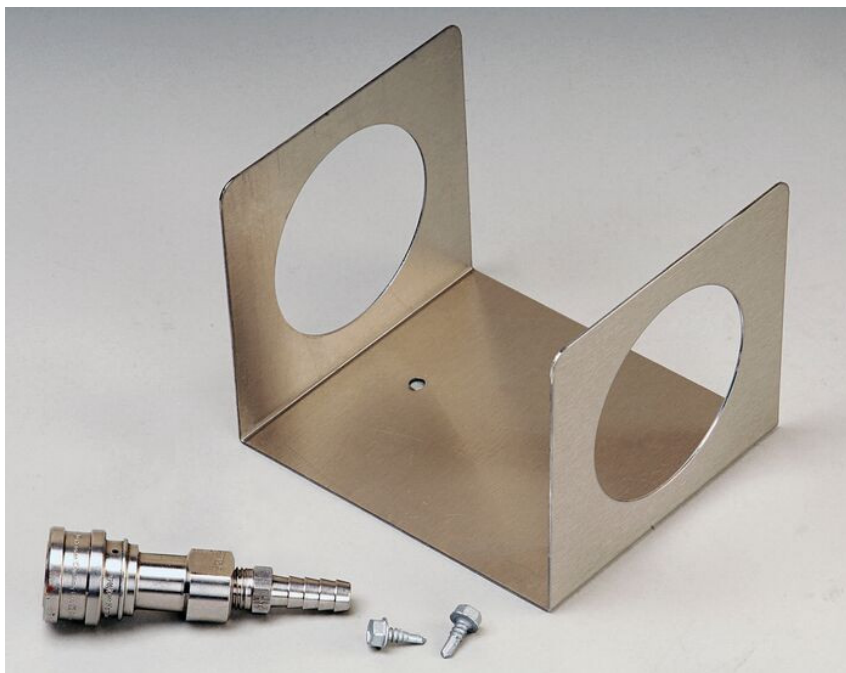

Short Description

- Ideal for applications using the Coalescing/Carbon Filter
- Prevents vapor release when replacing filters

Description

The Remote Wall Bracket Kit is ideal for applications using the Coalescing/Carbon Filter. Under typical conditions, these HPLC filters can be attached directly to the vapor release fitting on a HPLC Disposal Can. Where space is at a premium, it may be necessary to locate the filter in a remote location near the Disposal Can, such as a back wall. This accessory Bracket Kit is ideal for such circumstances. It consists of a 5" tall No. 304 stainless steel holder with a female quick-disconnect fitting and two self-tapping mounting screws. The disconnect fitting has a valve to provide a safe, fully closed system preventing vapor release when replacing filters. Fittings accommodate 3/8" OD (10mm) tubing, not included.

Specifications

Model No	28151
UPC	697841119437
International Model No.	28151
Color	Silver
Material Specifications	Stainless Steel
Dimensions, Exterior	5" H x 4" W
	127mm H x 102mm W
Net Weight, lbs	1.0
Net Weight, kg	0.5