
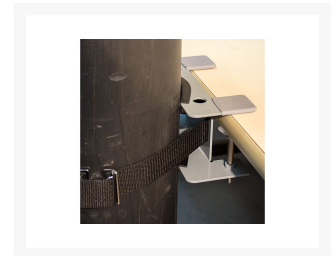


Steel Gas Cylinder Support Bracket, One 12" Cylinder Capacity, Bench Mount, Steel - 35258

Product Images

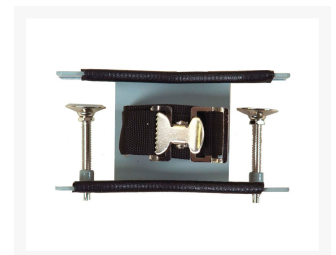

Short Description

- Brackets accept cylinder sizes ranging from 4- to 12-in (102- to 305-mm)
- Exclusive bumper guard edging with steel reinforced vinyl protects your cylinders from damage while providing extra grip
- Designed for bench mounting applications

Description

Cylinder support brackets provide a quick and easy solution for storing gas cylinders when space is at premium. Constructed from heavy gauge, powder-coated steel to ensure long service life and chemical resistance. Using dual 0.375-in (10-mm) diameter clamping posts with 1.5-in (38-mm) swivel pad faces, these models clamp to tables 1- to 2.75-in (26- to 70-mm) thick. Heavy-duty polypropylene support straps, 1.5- x 54-

in (38-mm x 1372-mm) and rated at 1138-lb (516-kg), feature steel cinch buckle and hold individual cylinders securely against bracket. For an improved level of safety for any installation, all brackets include openings to retrofit secondary support chains, one set for each cylinder station.

Helps comply with OSHA, NFP1, NFPA 55, and CGA.

Specifications

Model No	35258
UPC	697841210851
Cylinder Capacity	1
Minimum Cylinder Diam (in)	4
Cylinder Capacity with Max Cylinder Dia	(1) 12 inch
Mount Type	Bench
Strap Style	Poly Cinch Strap with Steel Buckle
Straps Per Cylinder	1
Material Specifications	Steel
Complies with NFPA 1 Fire-Code	Yes
NFPA 55	Yes
OSHA Compliance	Yes
CGA	Yes
Net Weight, lbs	3
Net Weight, kg	1.4
Net Dimensions (W x D x H)	8 x 4 x 4.25 inches
	203 x 102 x 108 mm