


Flammable Vertical Drum Cabinet

Product Images





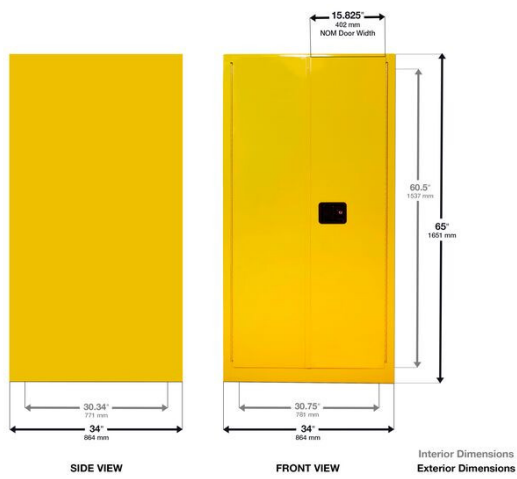

Short Description

- **Sturdy Construction:** Welded, 18-gauge steel with double-wall construction and 1.5 inches of insulating air space provides excellent fire resistance
- **Spill Containment:** Reinforced liquid-tight sump contains spills, ensuring safety and compliance
- **Enhanced Safety:** Dual vents with flame arresters provide additional safety features
- **Convenient Storage:** Removable, adjustable half-depth shelf accommodates drums with pumps or funnels
- **Drum Handling:** Certain models include drum rollers to facilitate easy movement of heavy drums
- **Configuration Options:** Available in multiple colors and various drum capacities to meet specific needs
- **Compliance:** Fully compliant with OSHA and NFPA regulations

Description

Eagle Vertical Drum Safety Cabinets are designed to safely store larger quantities of hazardous materials, such as flammable solvents, in compliance with OSHA and NFPA regulations. Constructed from sturdy 18-gauge steel with double-wall construction and 1.5 inches of insulating air space, these cabinets provide excellent fire resistance. They feature a reinforced liquid-tight sump to contain spills and dual vents with flame arresters for enhanced safety. A removable, adjustable half-depth shelf allows for additional storage and accommodates drums with pumps or funnels, making them suitable for use as pumping or waste accumulation stations. Certain models also include drum rollers, which facilitate the movement of heavy drums, making it easier to load and unload them from the cabinet.

To configure your safety cabinet, select the appropriate drum capacity, ranging from one drum to two 55- or 60-gallon drums, to fit your storage requirements. Choose from different colors to designate the type of materials stored, such as yellow for flammables, red for combustibles, gray or white for waste, or brown/neutral for labs. Additionally, choose between manual closing or self-closing doors to ensure

compliance with your safety standards. Self-closing doors offer enhanced safety by automatically latching in emergency situations, minimizing human error and ensuring compliance.

Eagle Vertical Drum Safety Cabinets are designed to meet the complex requirements of handling and storing hazardous materials, providing a reliable solution for industries such as laboratories, manufacturing facilities, and chemical processing plants. By configuring your cabinet, you ensure a safe and compliant storage solution tailored to your unique needs.

Specifications

Model No	DRV
----------	-----

Prop 65

⚠️WARNING: This product can expose you to chemicals including Lead, which is known to the State of California to cause cancer and birth defects or other reproductive harm. For more information go to www.P65Warnings.ca.gov.

Product Options

Color:	Red
	Yellow
	Gray
	Beige
	White
Capacity:	55 gallon
	110 gallon
Door Style:	Self-close
	Manual Close
	Bi-Fold